


β -Blockers High pH


Columns:	Listed on chromatograms
Dimensions:	150 × 4.6mm ID, 5 μm(all columns)
Mobile Phase:	MeOH: 5mM NH₄HCO₃(pH=10.0)=70:30
Flow Rate:	1.0ml/min
Temperature:	Ambient
Detection:	UV 220nm
Sample:	1. Pindolol 2. Metoprolol 3. Bisoprolol 4. Propranolol 5. Alprenolol

*ZORBAX,HC,TC,Eclipse,Hypersil are registered trademarks of Agilent Technologies.Gemini,Luna are registered trademarks of Phenomenex.Symmetry,XBridge are registered trademarks of Waters Corporation.Ineretsil is registered trademark of GL Sciences Inc.Kromasil is registered trademark of Eka Chemicals.BECKMAN Coulter is registered trademark of Beckman Coulter, Inc. SMT is registered trademark of Analtech Inc.Dikma Technologies is not affiliated with the above companies.


Spursil C18

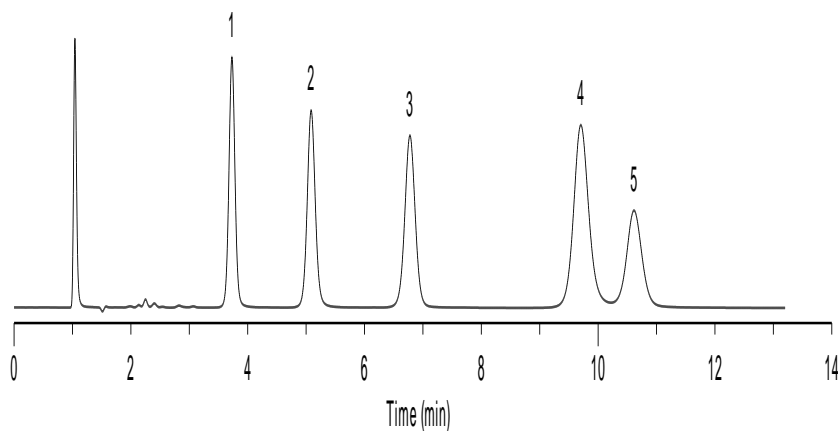


Spursil C18-EP

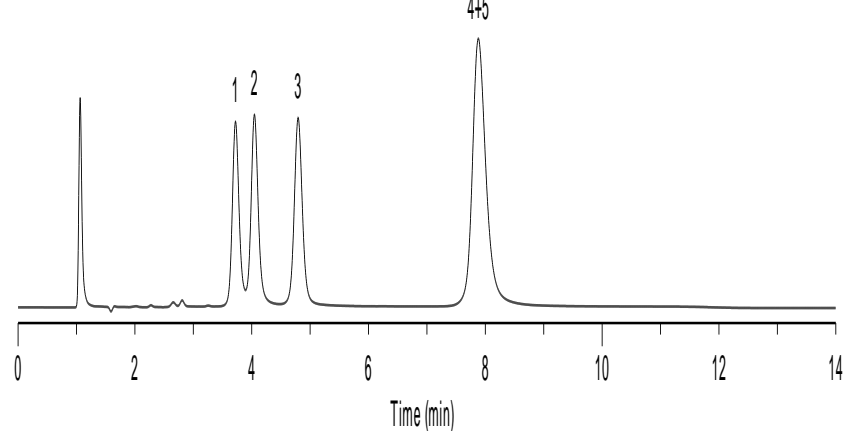


***Agilent ZORBAX Bonus-RP**

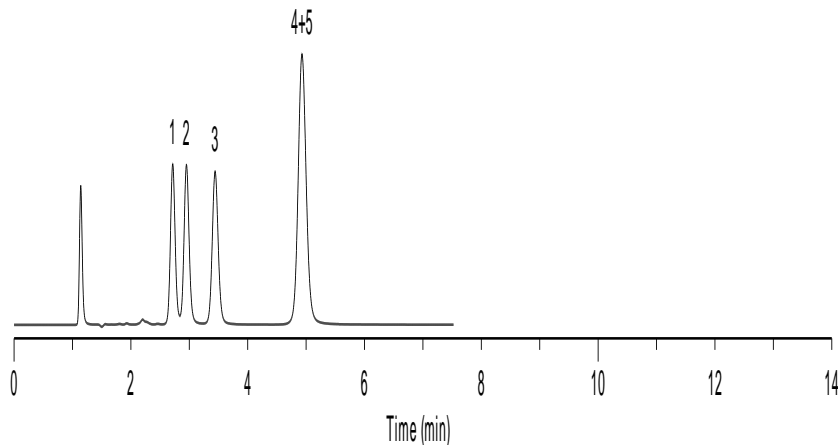


***Waters SymmetryShield RP18**

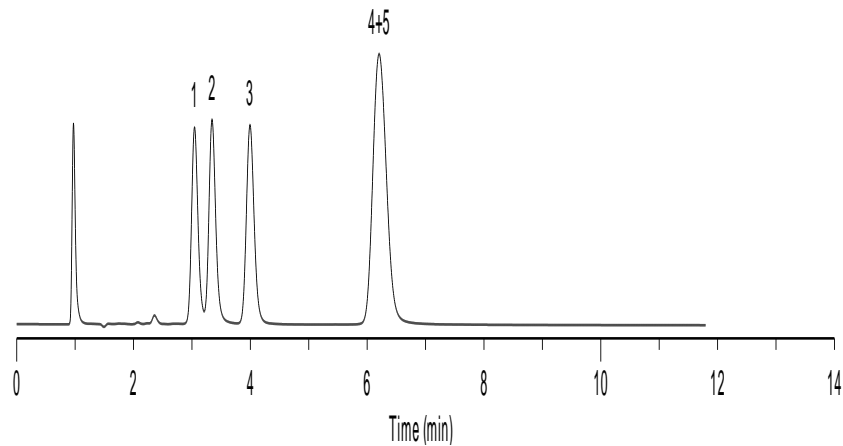


Diamonsil C18

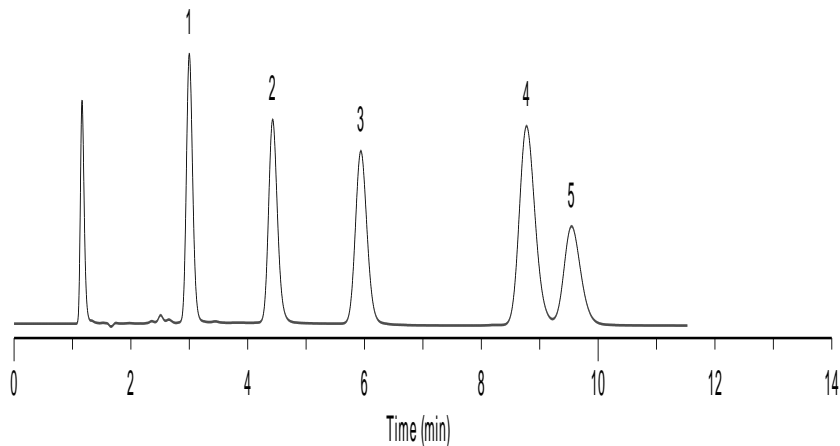


Inspire

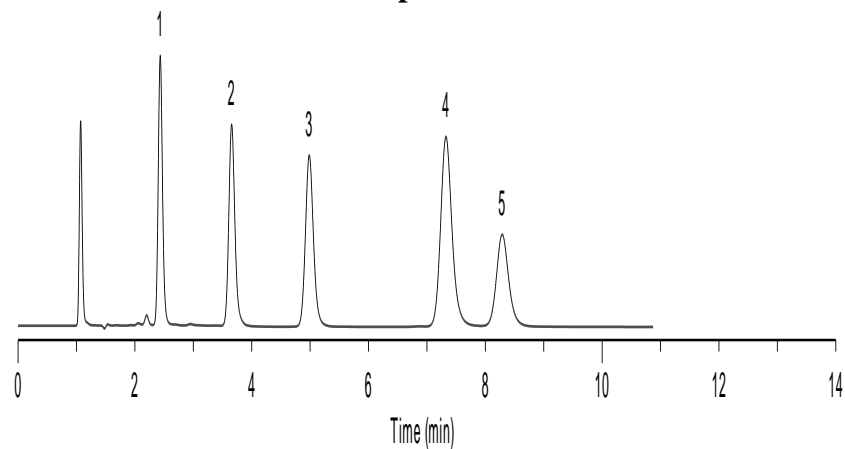


Platisil C18

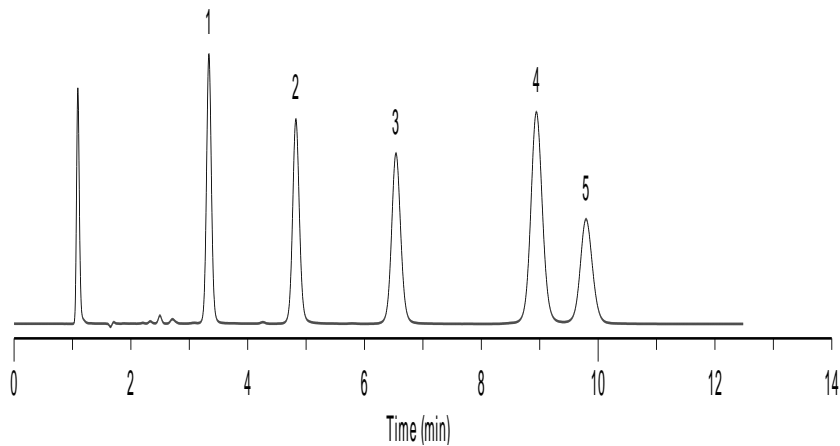


***Inertsil ODS-3**

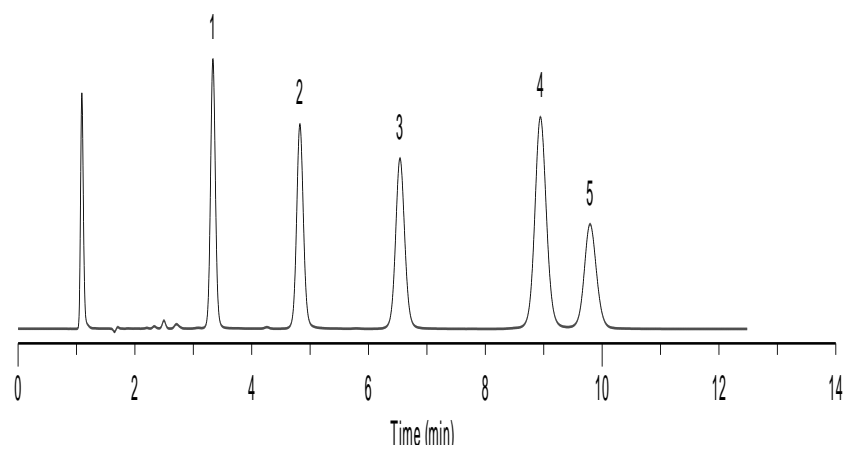


***Agilent HC-C18**

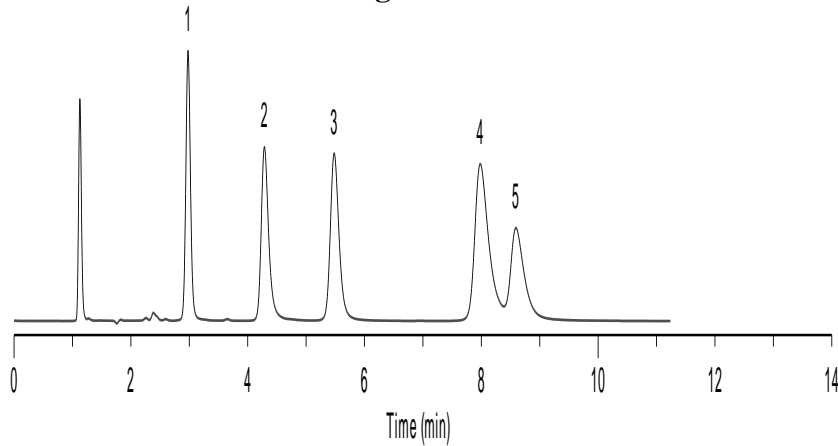


***Agilent TC-C18**

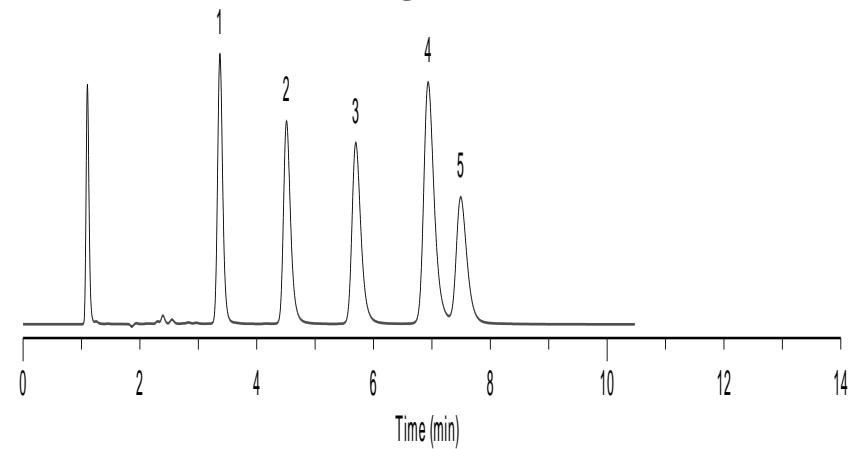


***Agilent ZORBAX SB-C18**

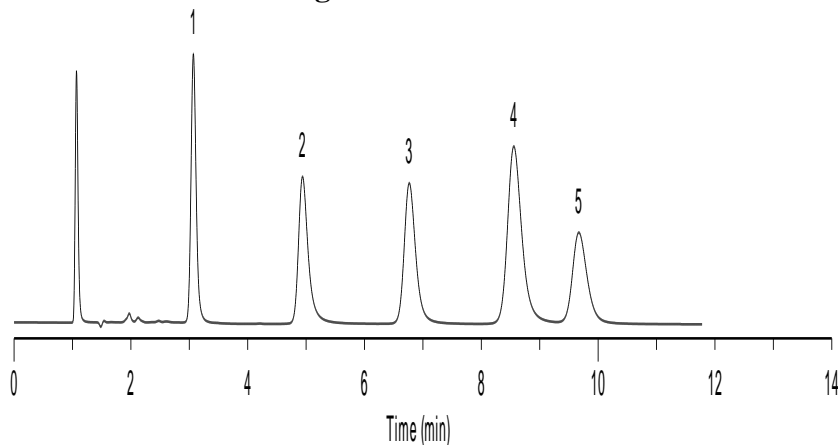


***Agilent Eclipse XDB-C18**

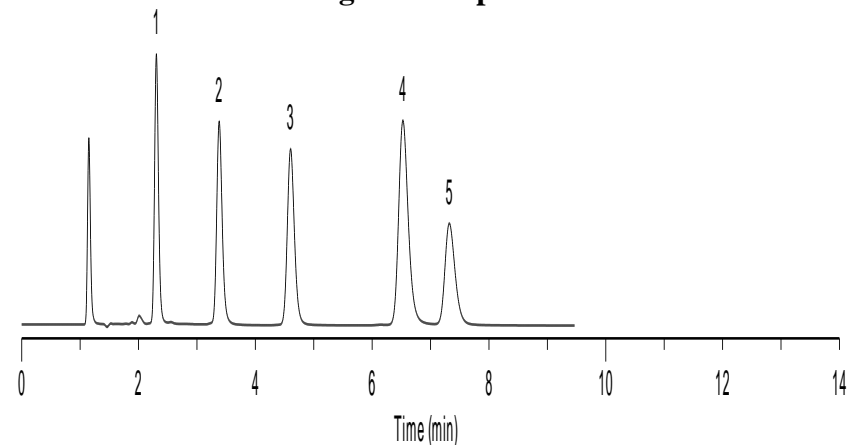


***Phenomenex Gemini C18**

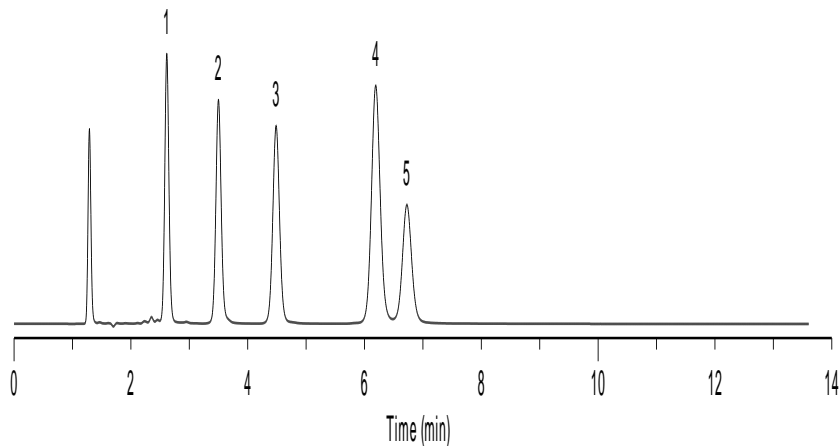


***Phenomenex Luna C18(2)**

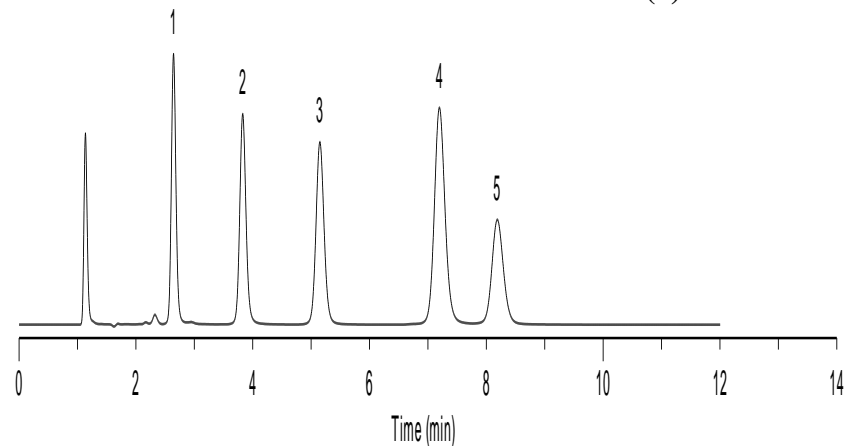


***Waters XBridge C18**

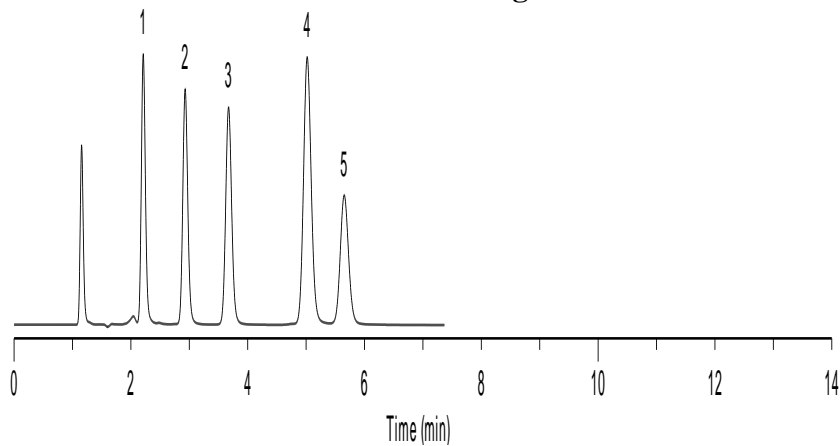


***Waters Symmetry C18**

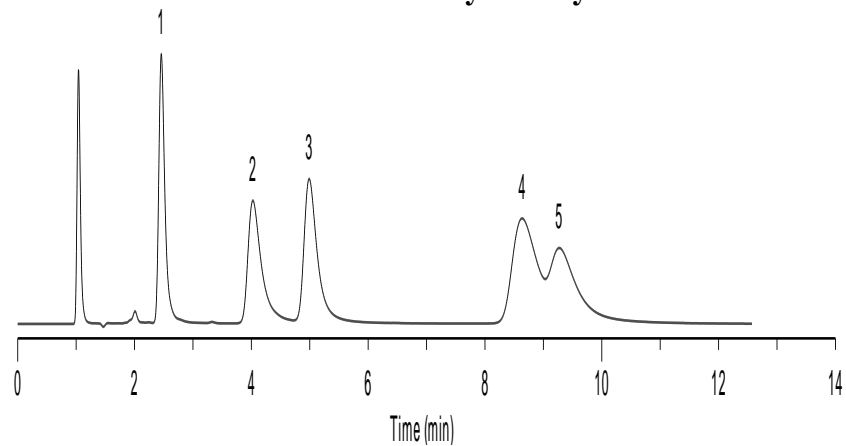


***Kromasil C18**

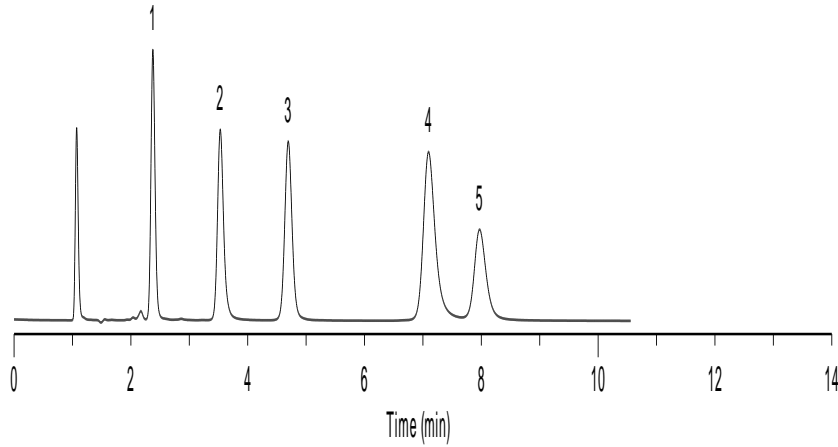


***Hypersil BDS C18**

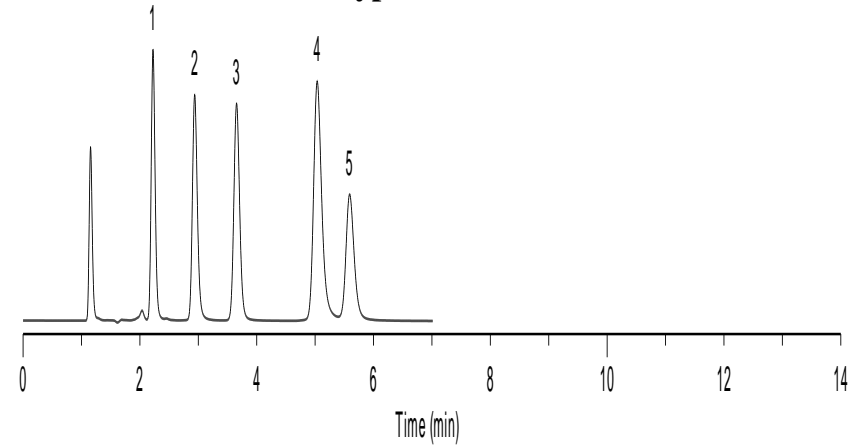

***BECKMAN Coulter C18**

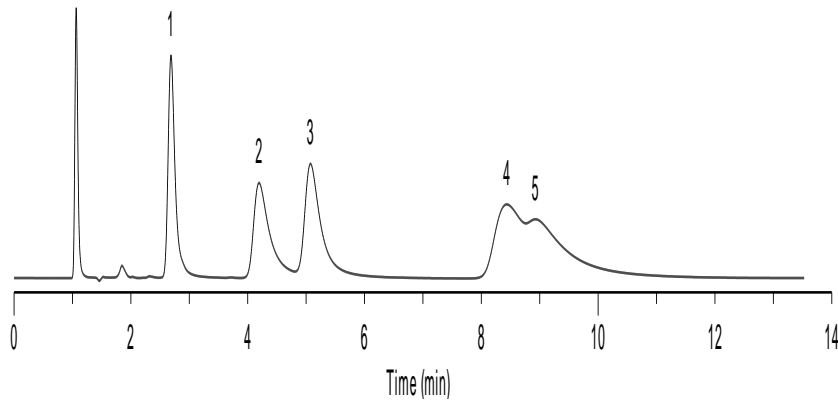

***SMT C18**

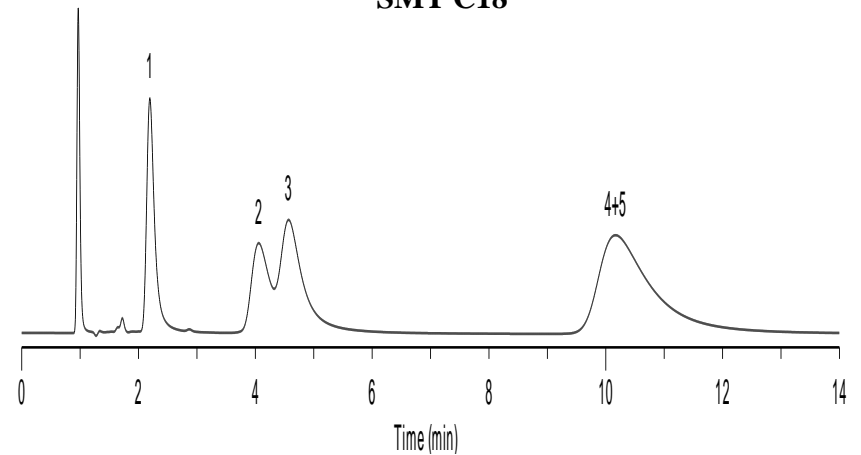

[China](#)

[Canada](#)

[US](#)