


Sample:  
Gel Filtration Standard  
(Produced by Bio-Rad)  
Injection volume = 20 $\mu$ l

1. Thyroglobulin
2.  $\gamma$  Globulin
3. Ovalbumin
4. Myoglobin
5. Cyanocobalamin

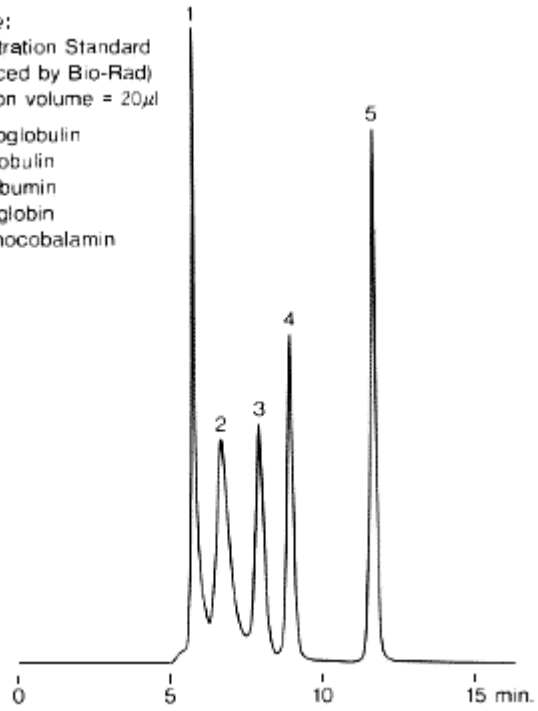

Column : Shodex PROTEIN KW-802.5 (8.0mmID\*300mm)  
Eluent : 50mM Phosphate buffer (pH7.0) + 0.3M NaCl  
Flow rate : 1.0mL/min  
Detector : UV (280nm)  
Column temp. : Room temp.