




Cadenza CD-C18

75 x 2 mm

Application

α - Solanine and α - Chaconine from Potato

ジャガイモ中の α -ソラニンと α -チャコニン

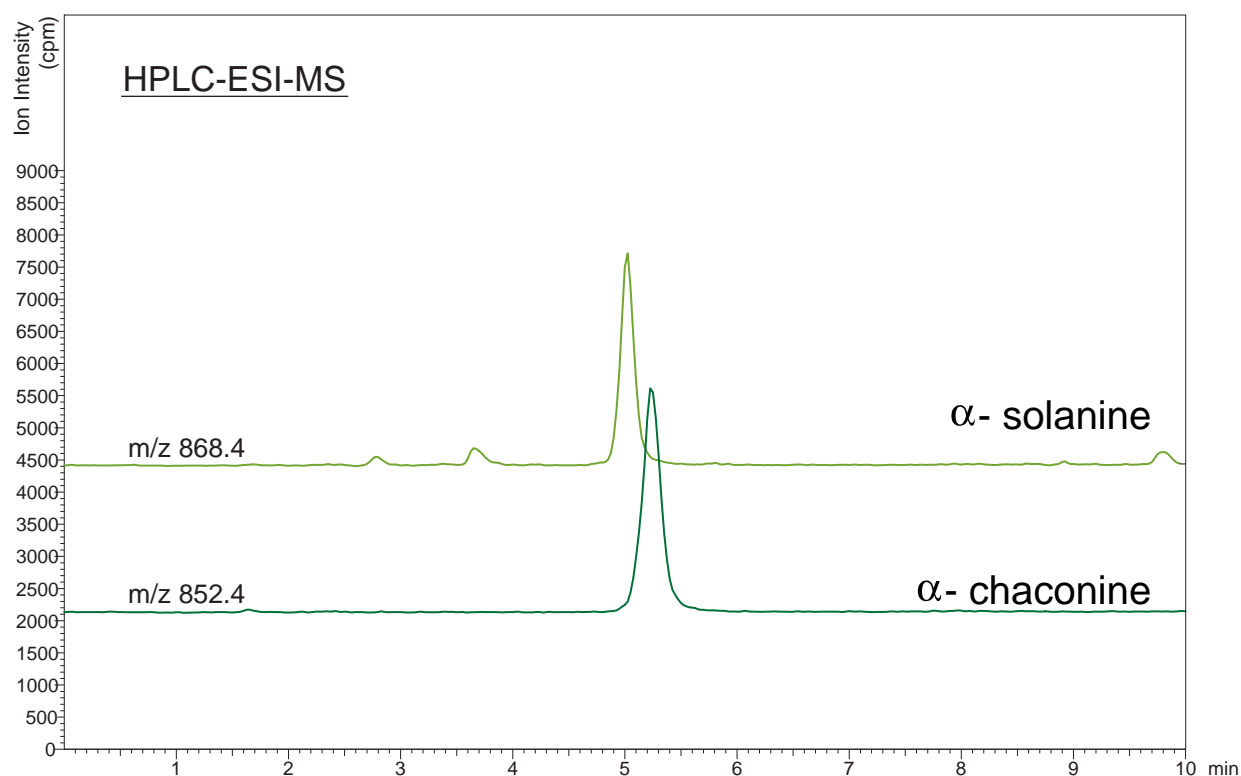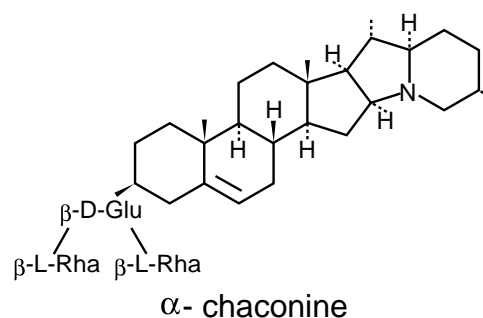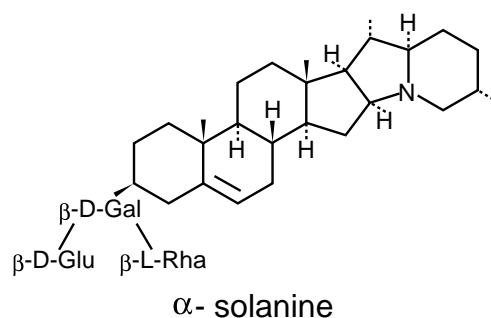

Cadenza CD-C18, 75 x 2 mm
 0.1% trifluoroacetic acid / acetonitrile = 75 / 35
 0.2 mL/min, 35 deg.C,
 LCMS-2010A, ESI positive

Courtesy of Dr. Fumio MATSUDA and Prof. Hisashi MIYAGAWA,
 Department of Agriculture, Kyoto University.