


Unison UK-C8

50 x 2 mm



Application

LC-MS Application for Antiseptic Substance

LC-MSアプリケーション (消毒薬, クロルヘキシジン)

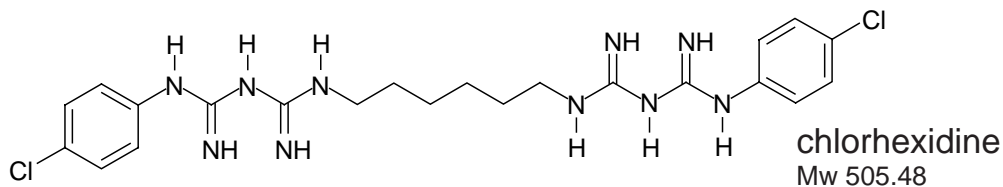

API4000 (TAKARA BIO INC.)

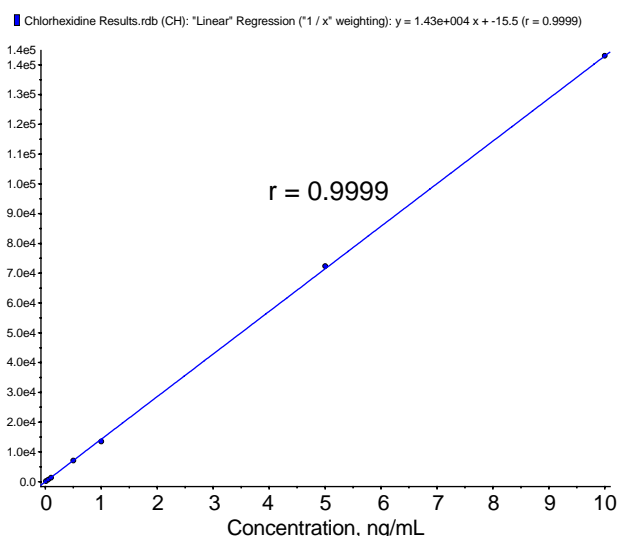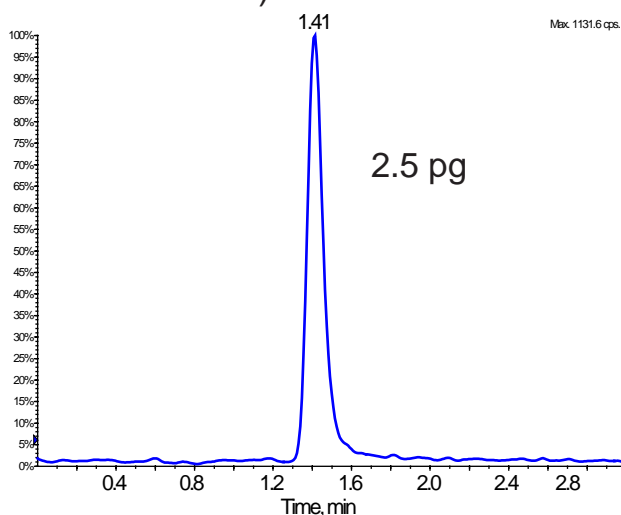

Unison UK-C8, 50 x 2 mm

0.1% formic acid / acetonitrile = 2 / 1, 0.2 mL/min, 5uL(500 fg/uL)

API4000: ESI, MRM Positive, Q1/Q3 : 505 / 336

Courtesy of J.Watanabe, TAKARA BIO INC.