





Scherzo SM-C18

250 x 4.6 mm

Application

Nojirimycin relative compounds

ノジリマイシン関連化合物

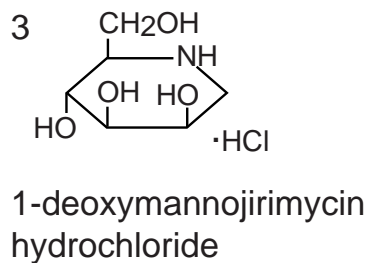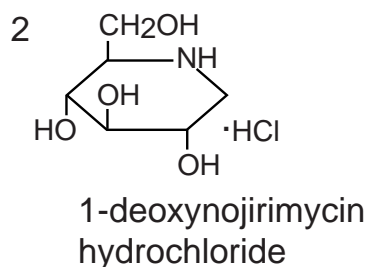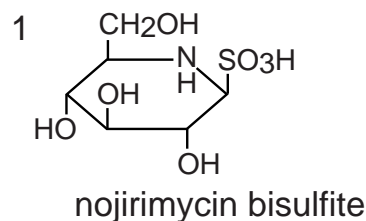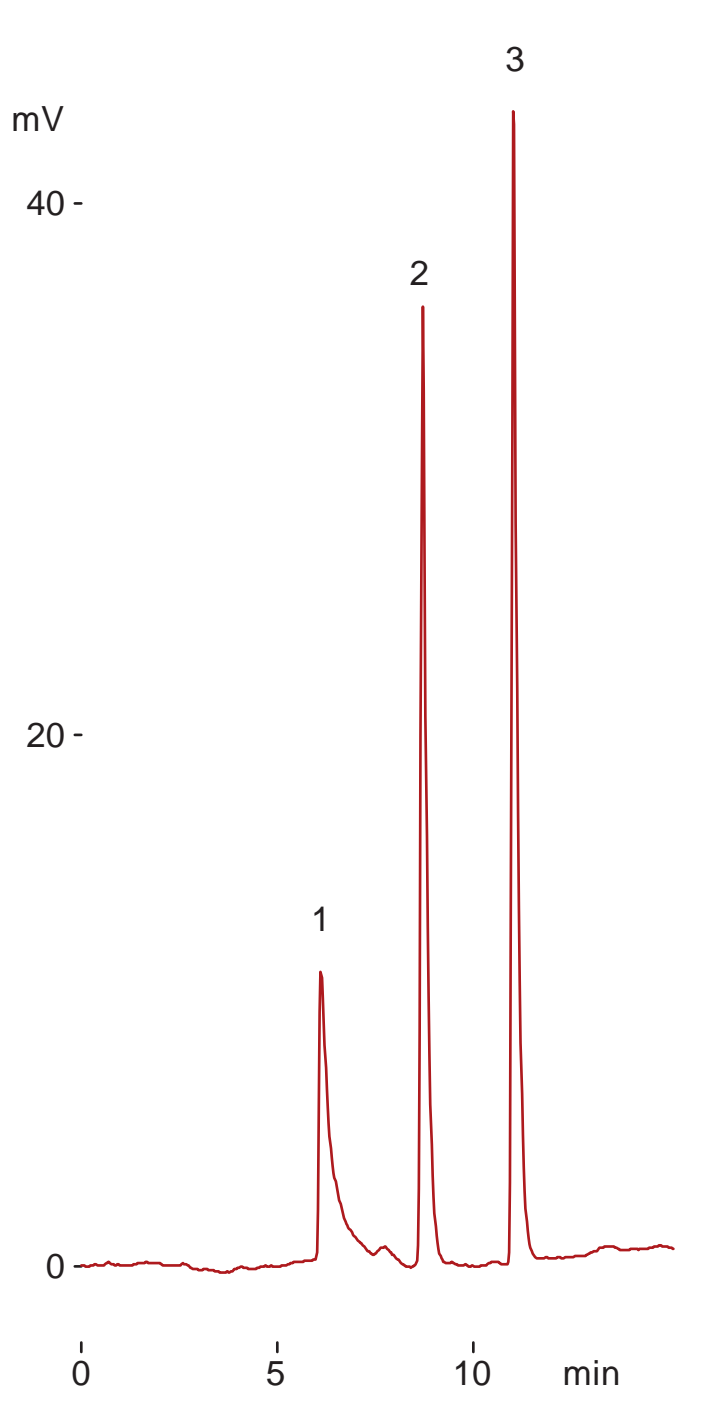

Scherzo SM-C18, 250 x 4.6 mm
 A: 5mM ammonium acetate
 B: 200mM ammonium acetate / acetonitrile = 80 / 20
 0 -40 %B (0-10min)
 0.5 mL/min (7MPa), 37 deg.C
 ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)
 5 uL (1.25-2.5ug)