

Scherzo SM-C18












150 x 4.6 mm

Application

Coumarins

クマリン類

A: water /acetic acid
= 100 /0.1
B: acetonitrile /acetic acid
= 100 /0.1

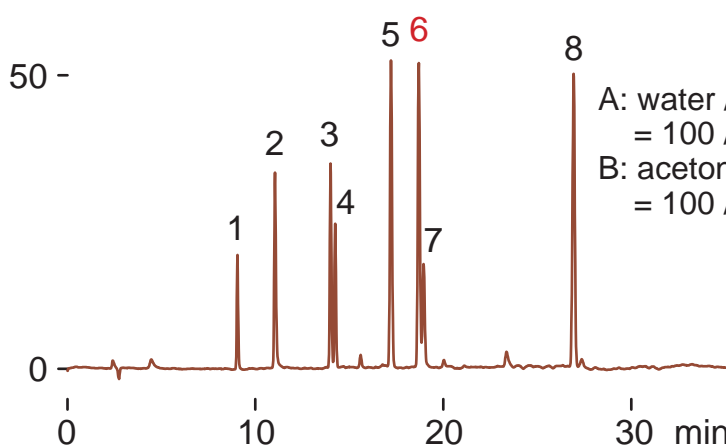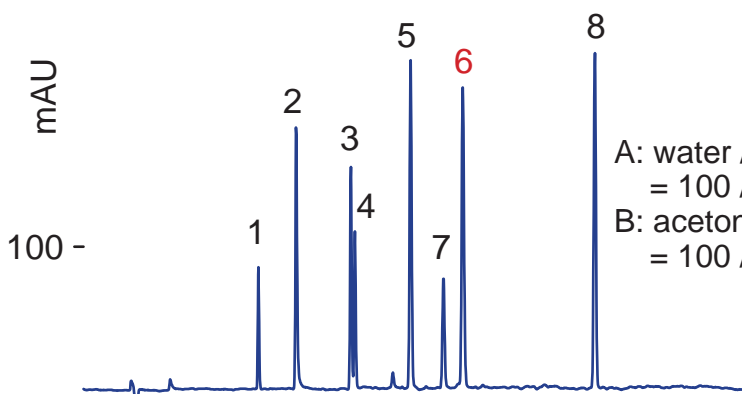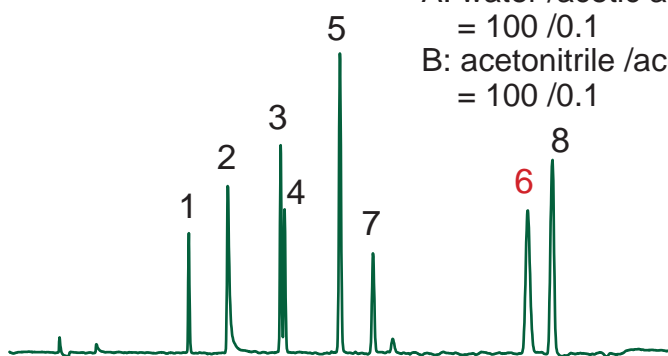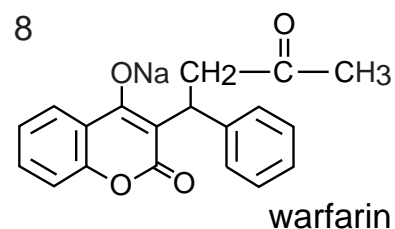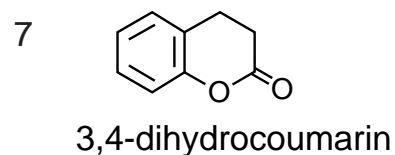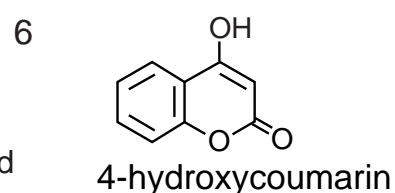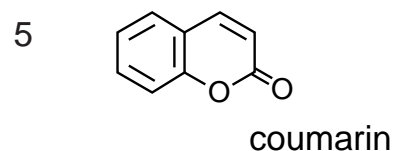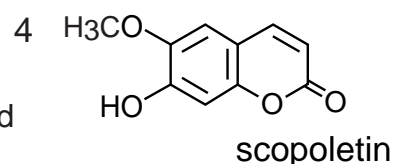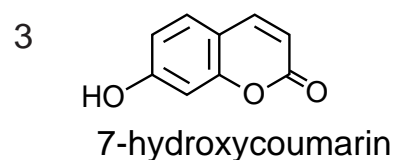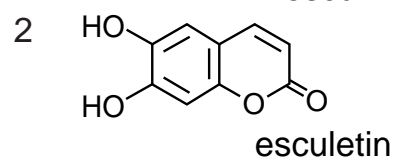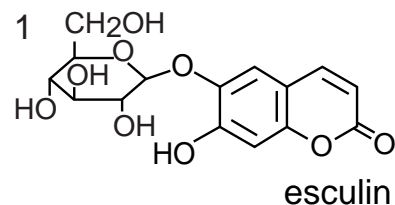

A: water /acetic acid
= 100 /0.5
B: acetonitrile /acetic acid
= 100 /0.5

A: water /acetic acid
= 100 /1
B: acetonitrile /acetic acid
= 100 /1

Scherzo SM-C18, 150 x 4.6 mm
0-60 %B (0-30min)
1 mL/min (9MPa), 37 deg.C, 270 nm
2 uL (0.2-1.2 ug)