


Presto FF-C18

150 x 4.6 mm

Application

Hyaluronic acid, chondroitin sulfuric acid

ヒアルロン酸とコンドロイチン硫酸の同時分析

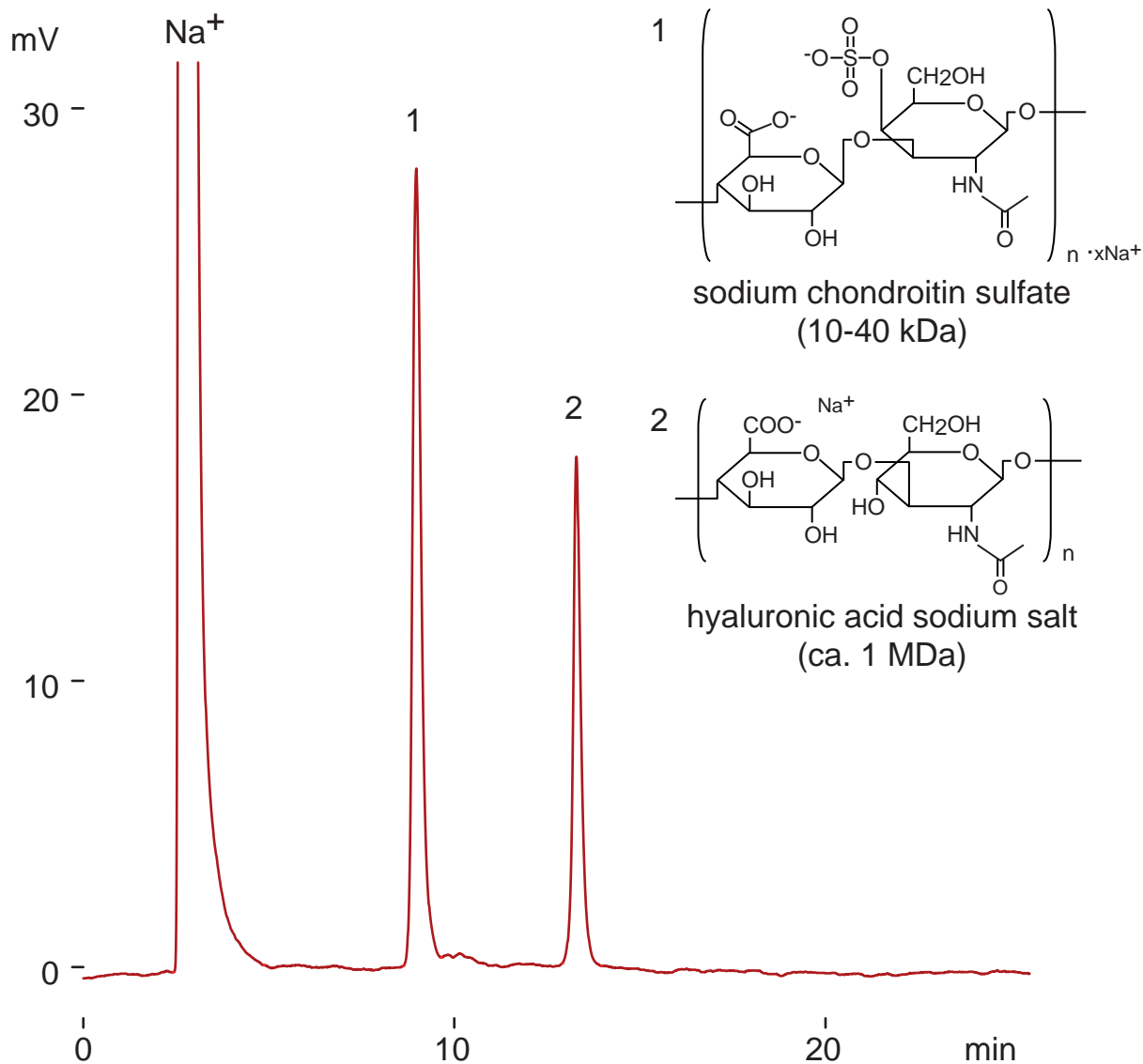

Presto FF-C18, 150 x 4.6 mm

A: water / trifluoroacetic acid = 100 / 0.05

B: acetonitrile / trifluoroacetic acid = 100 / 0.05

3-25 %B (0-15min), 3 %B (15-30min)

0.4 mL/min (20MPa), 25 deg.C

ELSD (spray chamber 10 deg.C, drift tube 50 deg.C)

20 uL (0.25% TFA-water)