


Milk Thistle Dry Extract using an UltraCore Column (II)

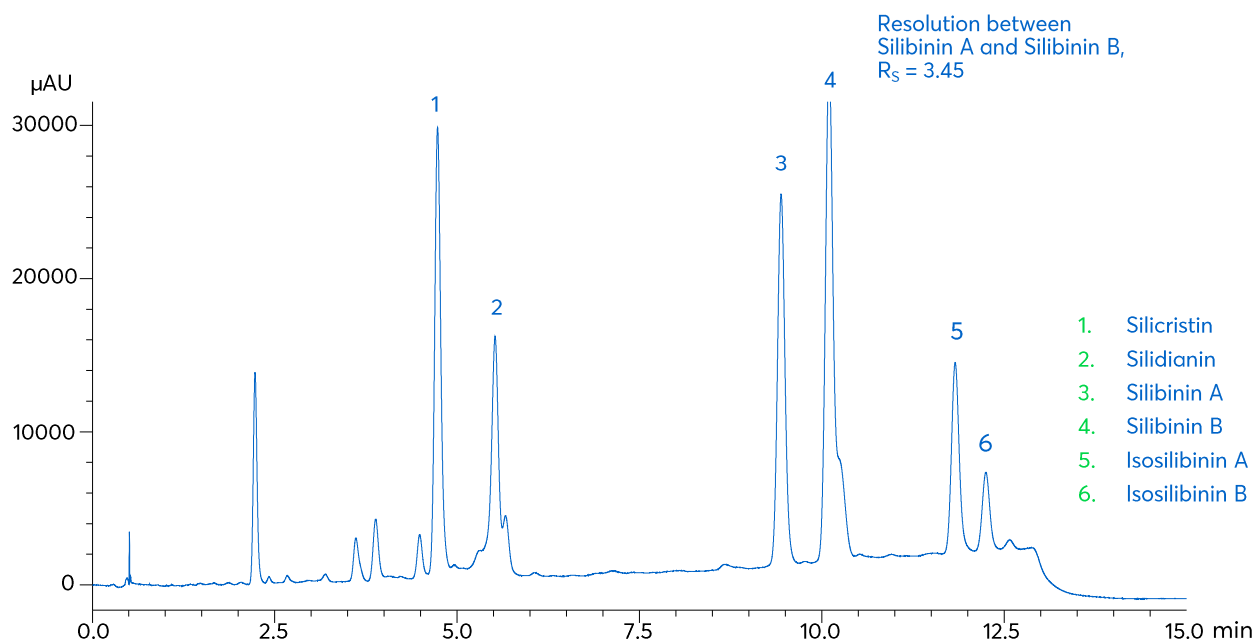

Method Details



CONDITIONS

Columns: Avantor® ACE® UltraCore 2.5 SuperC18
 Particle Size: 2.5 µm
 Dimensions: 75 x 4.6 mm
 Mobile Phase: A: MeOH/ H₂O/ Orthophosphoric acid (35/65/0.5)
 B: MeOH/ H₂O/ Orthophosphoric acid (50/50/0.5)

Gradient:

Time (mins)	%B
0.00	0
12.50	100
12.51	0
15.00	0

Flow Rate: 1.3 mL/min
 Injection Volume: 2 µL
 Detection: UV, 288 nm
 Sample: Milk thistle dry extract HRS

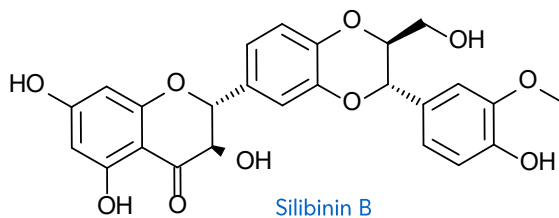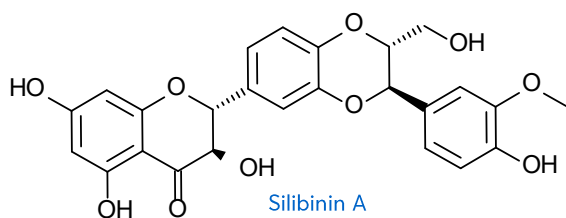

ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® UltraCore 2.5 SuperC18	HPLC Column	75 x 4.6 mm	CORE-25A-7546U

Reproduced with permission of Stephan Avramov, PhD, Chromassist Ltd, Bulgaria