
Determination of Honokiol and Magnolol in Huoxiangzhengqi Water

ACE®
Ultra-inert
UHPLC & HPLC Columns


Application #AN6180

Conditions


Column:	ACE Excel 5 C18
Dimensions:	250 x 4.6 mm
Part Number:	EXL-121-2546U
Mobile phase:	Methanol:Acetonitrile:Water (40:20:40 v/v/v)
Flow Rate:	1.0 mL/min
Injection:	10 µL
Temperature:	Ambient
Detection:	UV, 294 nm
Monograph method:	ChP. 2015, Vol. 1

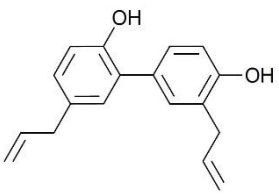

1. Honokiol

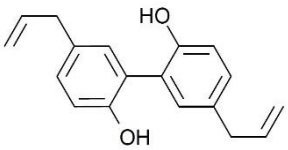

2. Magnolol

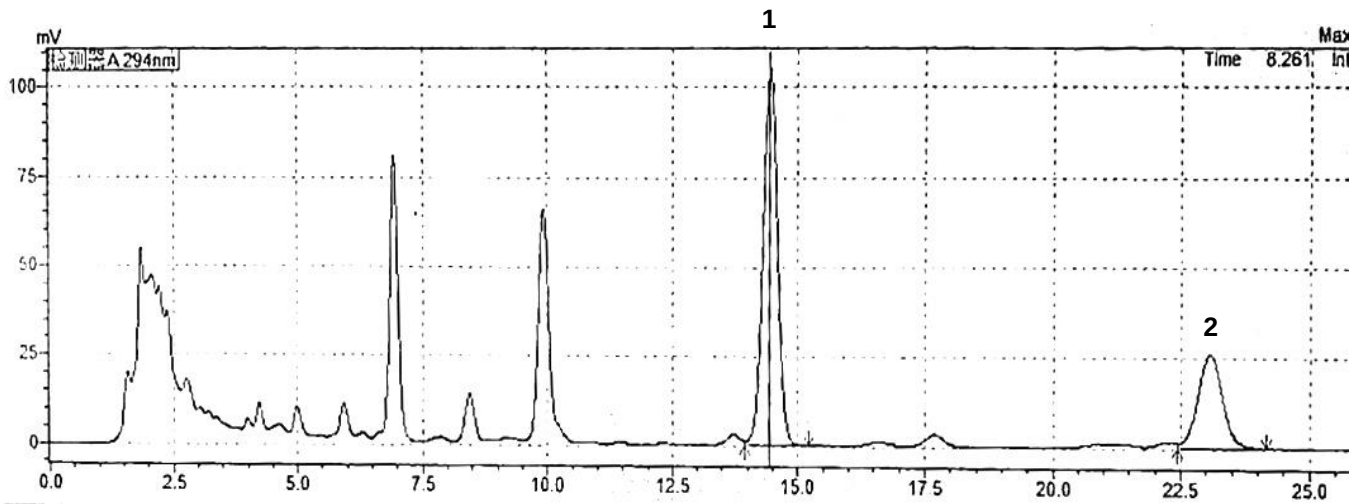

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.