


EP Assay for Milk Thistle Dry Extract

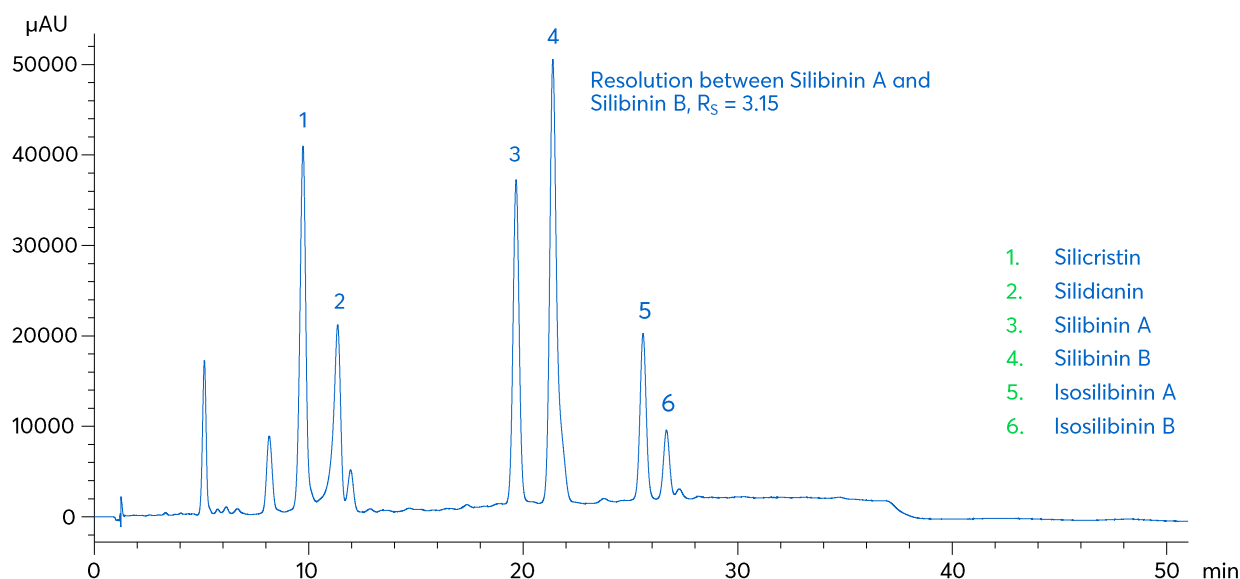

Method Details



CONDITIONS

Columns: Avantor® ACE® 5 C18
 Particle Size: 5 µm
 Dimensions: 125 x 4.0 mm
 Mobile Phase: A: MeOH/H₂O/Orthophosphoric acid (35:65:0.5 v/v/v)
 B: MeOH/H₂O/Orthophosphoric acid (50:50:0.5 v/v/v)

Gradient:

Time (mins)	%B
0	0
28	100
35	100
36	0
51	0

Flow Rate: 0.8 mL/min
 Injection Volume: 5 µL
 Detection: UV, 288 nm
 Sample: Milk thistle dry extract HRS
 Method: EP 10.0

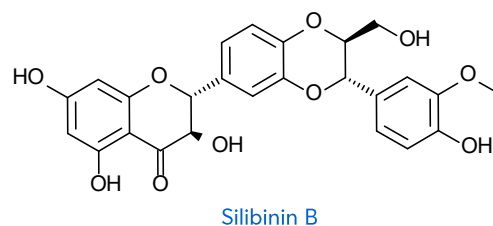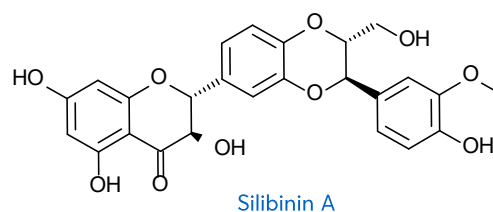

Always check the latest and official method information from the relevant pharmacopoeia prior to analysis

ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® 5 C18	HPLC Column	125 x 4.0 mm	ACE-121-1204

Reproduced with permission of Stephan Avramov, PhD, Chromassist Ltd, Bulgaria