

Chromatography Solutions


Application note #7780


Traditional Chinese Medicines Chemical Profiling

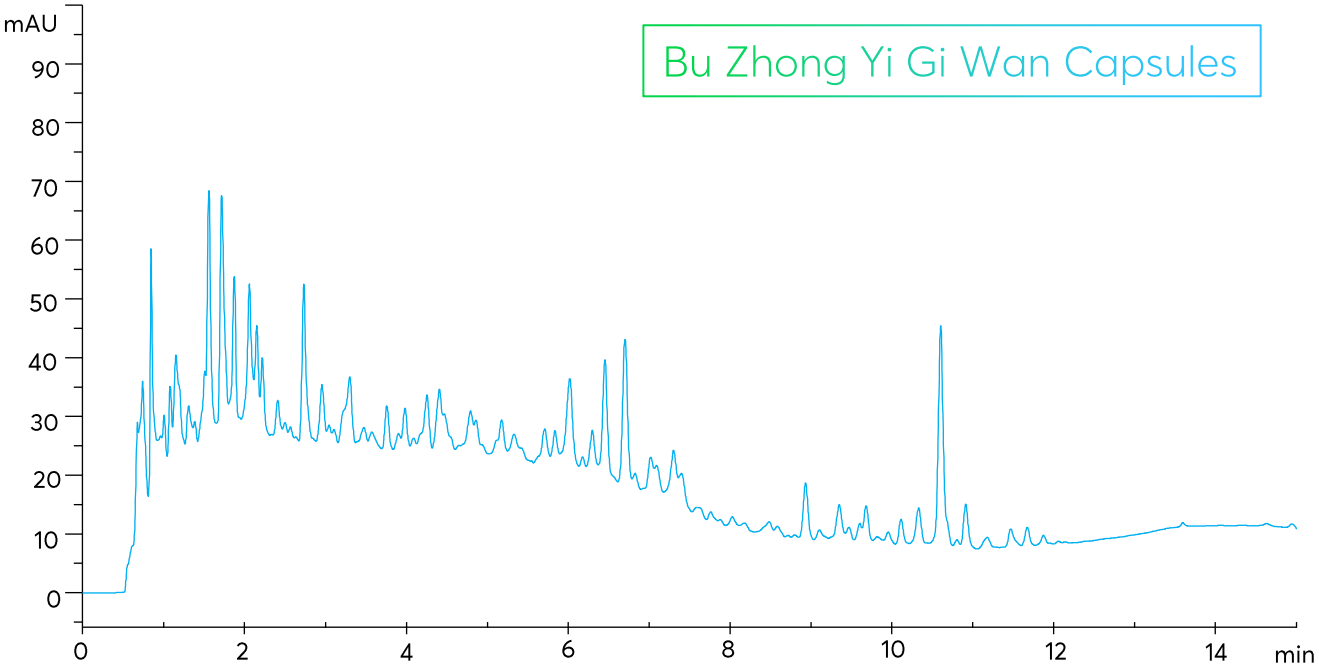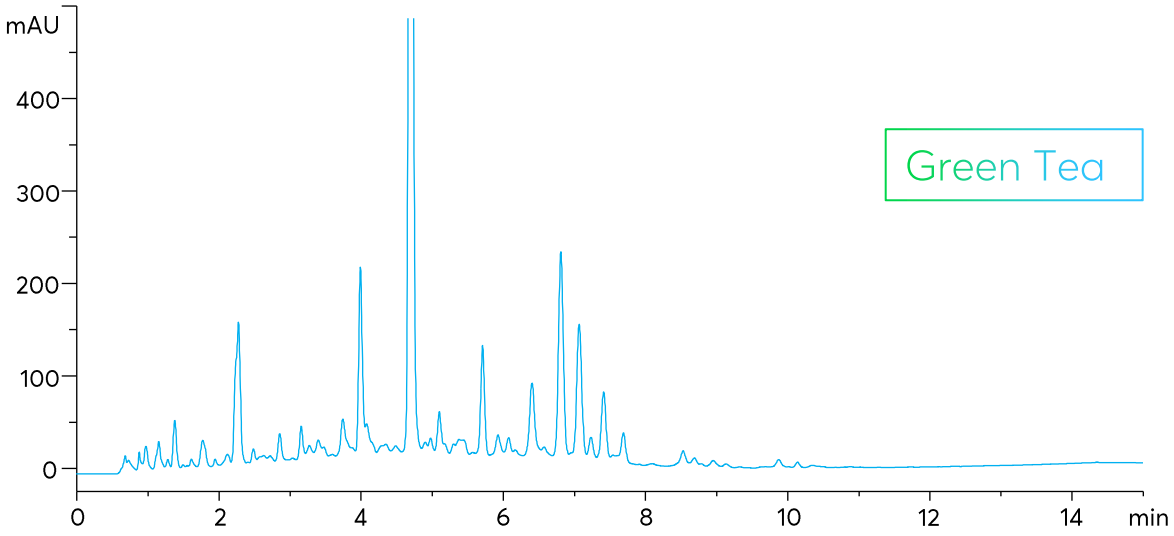

Avantor® ACE®


mAU
400
300
200
100
0
0
2
4
6
8
10
12
14
min
Green Tea

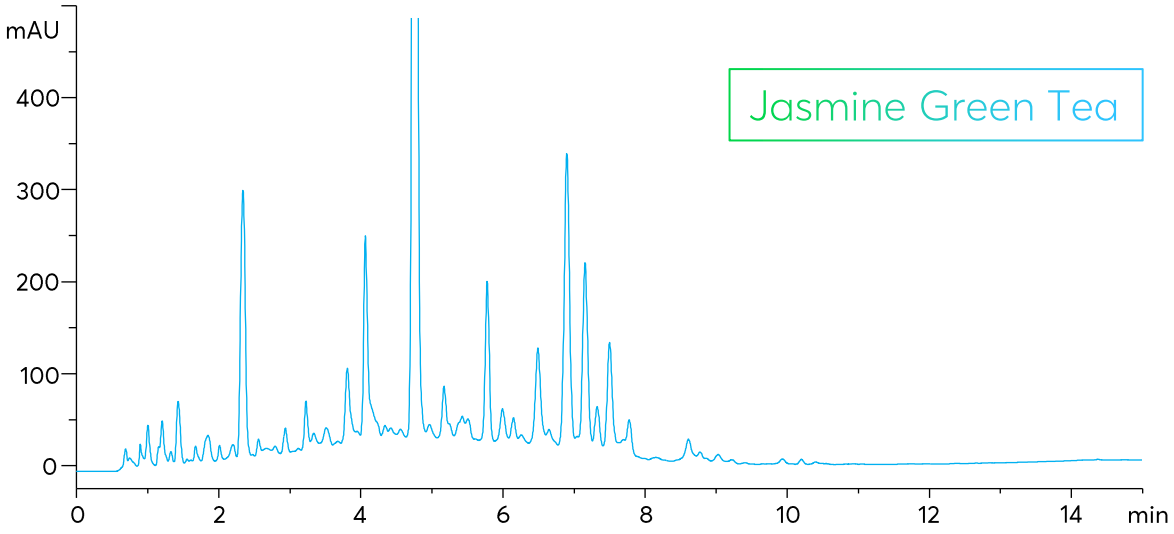
mAU
400
300
200
100
0
0
2
4
6
8
10
12
14
min
Jasmine Green Tea

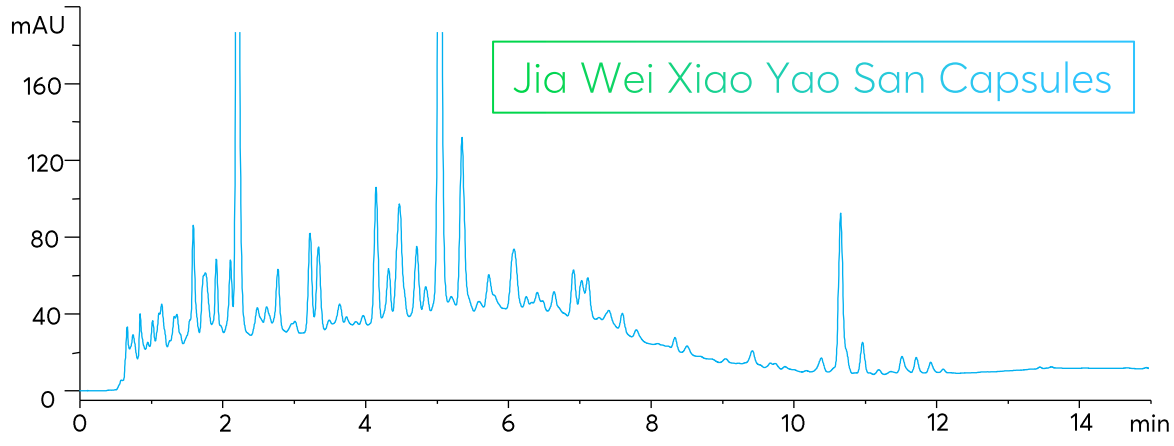
mAU
160
120
80
40
0
0
2
4
6
8
10
12
14
min
Jia Wei Xiao Yao San Capsules

Method Details

CONDITIONS

Columns: Avantor® ACE® Excel® 5 C18

Particle Size: 5 µm

Dimensions: 100 x 4.6 mm

Mobile Phase: A: 0.1% formic acid in H_2O

B: 0.1% formic acid in MeOH

Gradient:

Time (mins)	% B
0	5
19	100
20	100
22	5
25	5

Flow Rate: 1.5 mL/min

Injection: 10 µL

Temperature: Ambient 21 °C

Detection: UV, 254 nm

Sample preparation: Teas: steeped in 30 mL of 80 °C water for 10 min

Capsules: 15 tablets ultrasonicated in 15 mL of 80 °C water for 10 min

Samples were cooled to ambient and filtered via syringe filters before injection

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® Excel® 5 C18	HPLC Column	100 x 4.6 mm	EXL-121-1046U
VWR® Syringe Filters	Syringe filter nylon 25 mm 0.45 µm non-sterile	Pack of 100 items	514-1267
Methanol	VWR CHROMANORM® gradient grade	2.5 L	20864.320
Water	VWR CHROMANORM® super gradient grade	2.5 L	83650.320
Formic acid	VWR NORMAPUR® analytical reagent	100 mL	20318.180

Reproduced with permission of Prof. Andrew Shalliker, Chromaspeed Pty Ltd, and Western Sydney University, North Parramatta, Australia.

vwr.com/ace
chromsupport@avantorsciences.com