
Determination of Artesunate for Injection

ACE®
Ultra-inert
UHPLC & HPLC Columns


Application #AN6200

Conditions

Column:	ACE Excel 3 C18
Dimensions:	100 x 4.6 mm
Part Number:	EXL-111-1046U
Mobile phase:	10 mM potassium dihydrogen phosphate pH 3.0:MeCN (56:44 v/v)
Flow Rate:	1.0 mL/min
Injection:	20 µL
Temperature:	30 °C
Detection:	UV, 216 nm
Monograph method:	ChP. 2015

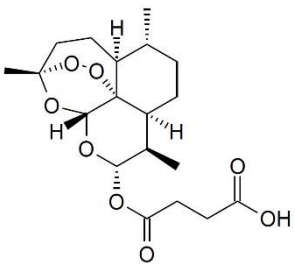

1. Artesunate

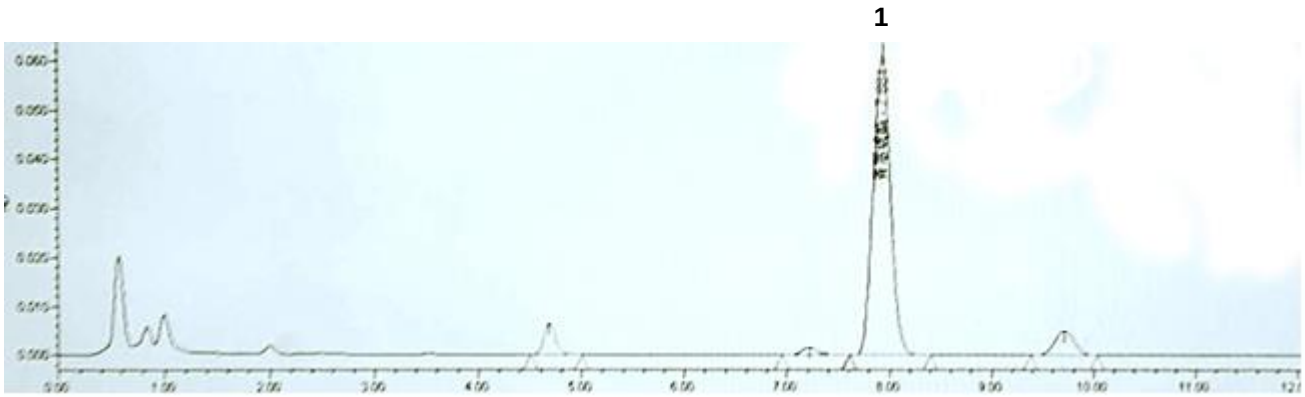

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.