


Isocratic analysis of atrazine herbicides

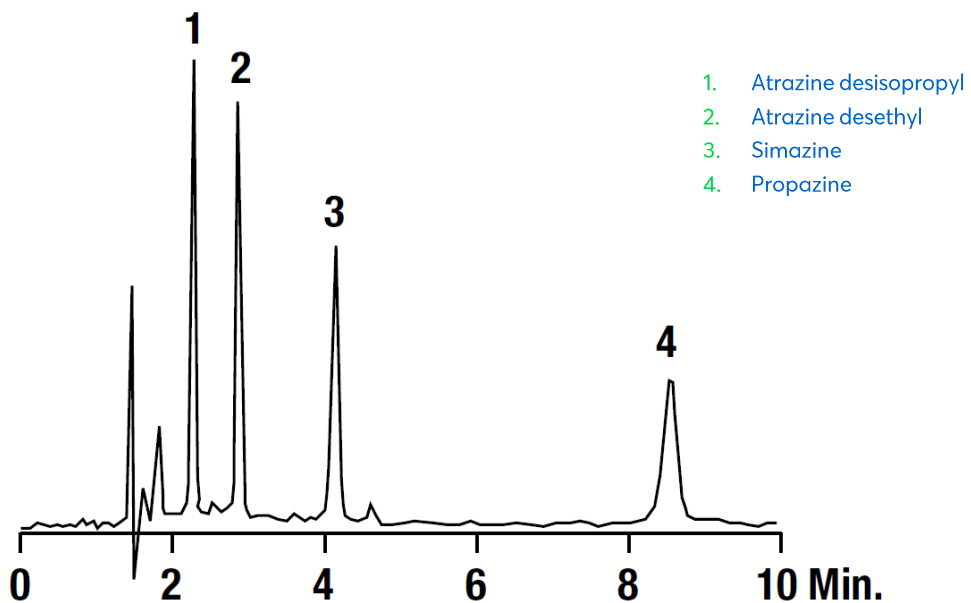

Method Details

CONDITIONS

Column:	Avantor® Genesis Cyano
Particle Size:	4 µm
Dimensions:	150 x 4.6 mm
Mobile Phase:	0.05% TFA in H ₂ O (pH3.0):MeCN
Flow Rate:	1.0 ml/min
Detection:	UV, 254 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® Genesis Cyano-RP	HPLC Column	100 x 4.6 mm	FM15965ER
Acetonitrile	VWR CHROMANORM® super gradient grade	2.5 L	83639.320
Water	VWR CHROMANORM® super gradient grade	2.5 L	83650.320
Trifluoroacetic acid	VWR HiPerSolv CHROMANORM® for LC/MS	10 x 1 mL	85049.001