


Artificial Sweeteners (Stevia Glycosides)

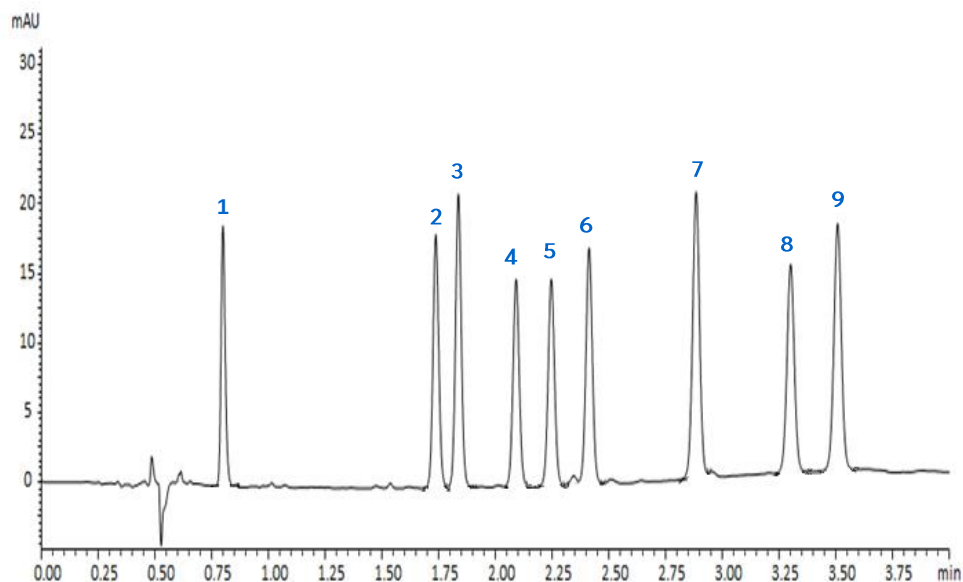

1. Rebaudioside D
2. Rebaudioside A
3. Stevioside
4. Rebaudioside F
5. Rebaudioside C
6. Dulcoside A
7. Rubusoside
8. Rebaudioside B
9. Steviolbioside

Method Details



CONDITIONS

Column: Avantor® ACE® Excel 2 SuperC18
 Particle Size: 2 µm
 Dimensions: 150 x 2.1 mm
 Mobile Phase: A: 10 mM sodium dihydrogen phosphate pH 2.8 in H₂O
 B: 10 mM sodium dihydrogen phosphate pH 2.8 in H₂O/MeCN (20:80 v/v)

Gradient:

Time (mins)	% B
0	39.5
4	48.0

Flow Rate: 0.6 mL/min
 Injection: 1 µL
 Temperature: 50 °C
 Detection: UV, 200 nm

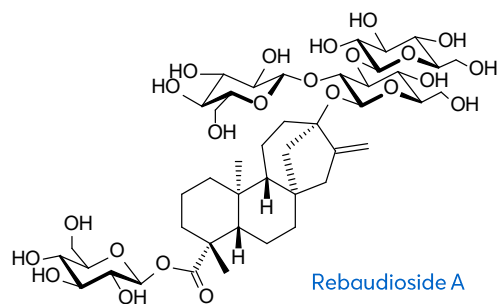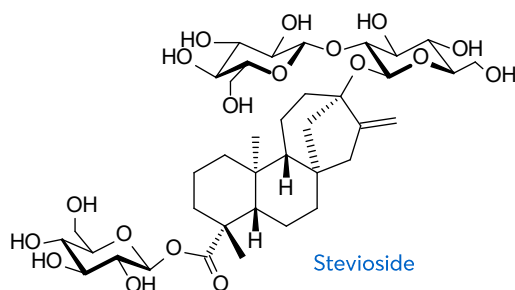

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® Excel 2 SuperC18	UHPLC Column	150 x 2.1 mm	EXL-1011-1502U

Reproduced with permission of Shimadzu Europe GmbH. Analysed on Shimadzu Nexera X2.
<http://images2.advanstar.com/PixelMags/lcgc-eu/pdf/2013-10-sp.pdf>