
Neomycin Sulfate and Sulfacetamide Sodium Separation


ACE®
Ultra-inert
UHPLC & HPLC Columns

Application #AN4510


Conditions

Column:	ACE 3 C18-PFP
Dimensions:	150 x 4.6 mm
Part Number:	ACE-1110-1546
Mobile Phase:	0.5% v/v TFA in H_2O/MeOH (97:3 v/v)
Flow Rate:	0.6 mL/min
Injection:	3 µL
Temperature:	30 °C
Detection:	ELSD, 30 °C
Instrument:	Chromaster

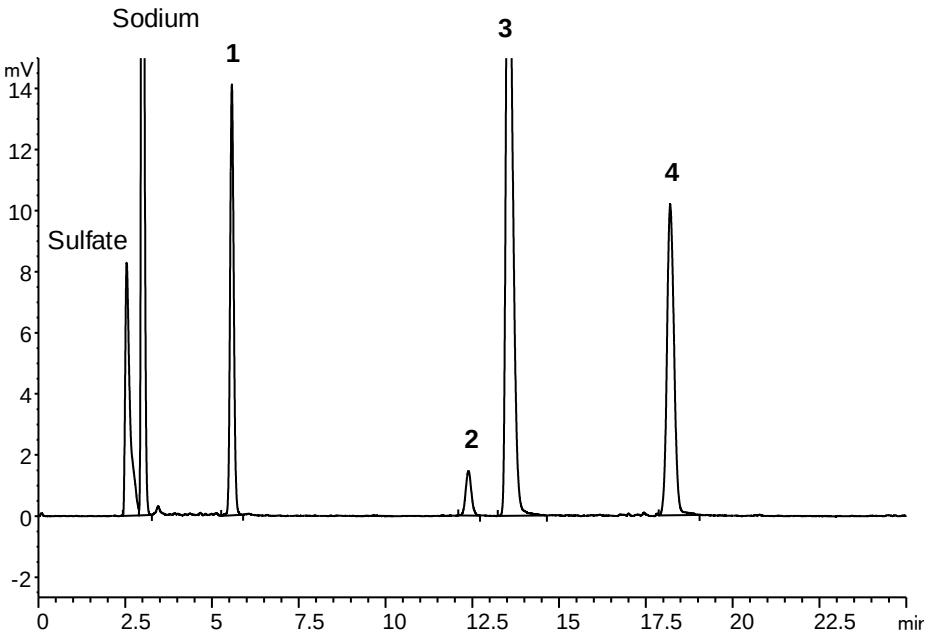


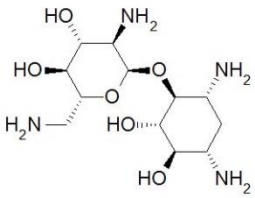


1. Neamine

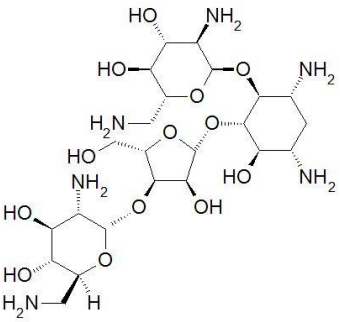

2. Neomycin C

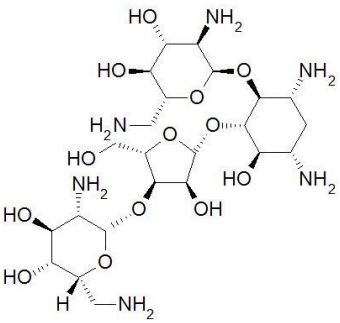

3. Neomycin B

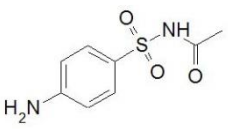

4. Sulfacetamide

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.