Amino Acids in Peas (*Pisum sativum*) by HPLC-HRAM-MS

ACE®
Ultra-inert
UHPLC & HPLC Columns


Application #**AN2660**

Conditions

Column:	ACE 3 AQ
Dimensions:	150 x 3.0 mm
Part Number:	ACE-116-1503
Mobile Phase:	A: 0.1% formic acid in H_2O
	B: 0.1% formic acid in MeCN

Time (mins)	%B
0	0
10	100

Flow Rate:	0.4 mL/min
Injection:	5 µL
Temperature:	30 °C
Detection:	Exactive Orbitrap high resolution MS
	ESI positive ion mode
	Capillary temperature: 350 °C

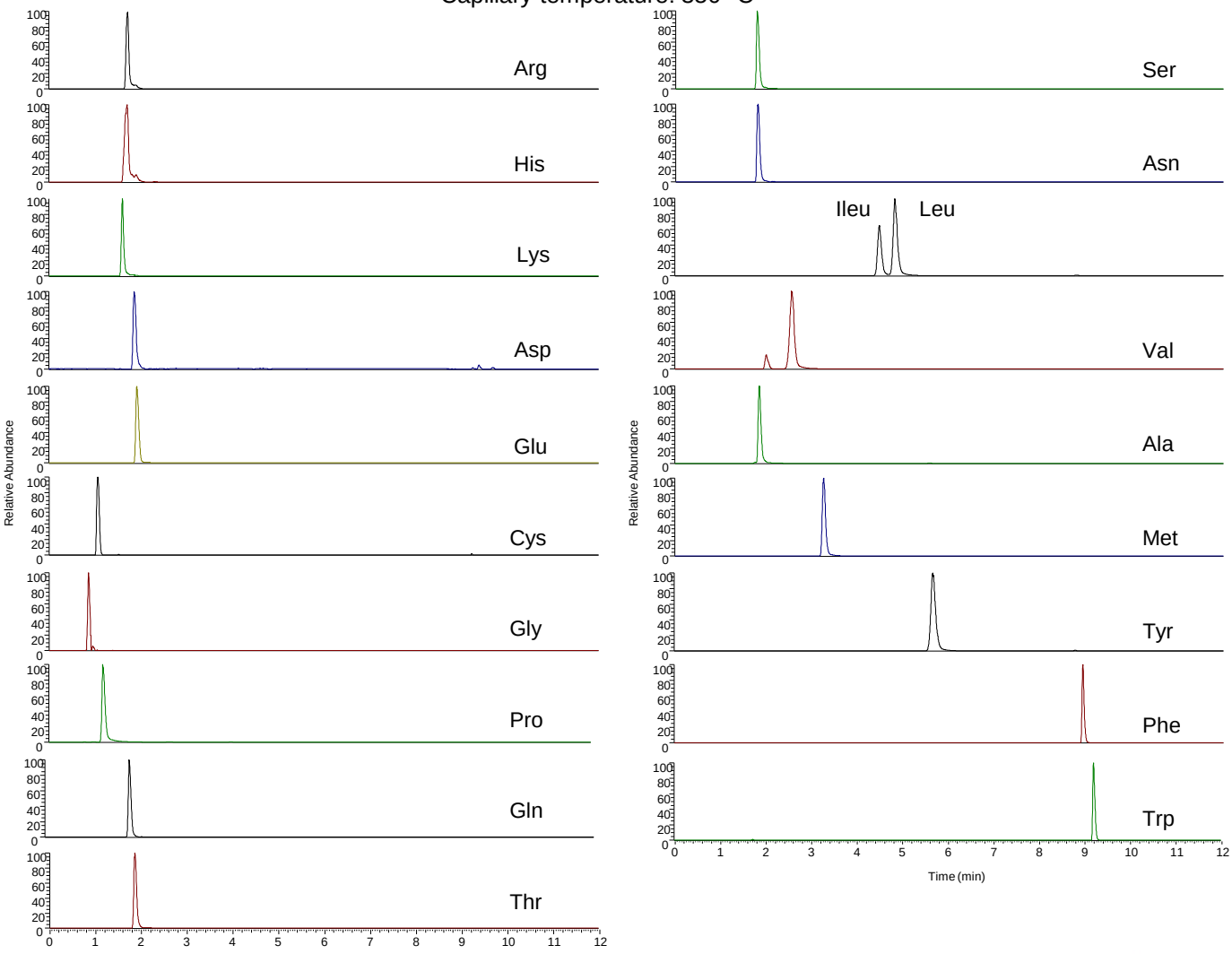

Reproduced with permission of Fera Science Ltd, York, UK
© 2015 ACT Ltd. All rights reserved. All products are available worldwide.