


1. Desmethyl venlafaxine
2. Hydroxybupropion
3. Hydroxyrisperidon
4. Bupropion
5. Zuclopenthixol
6. Fluphenazin

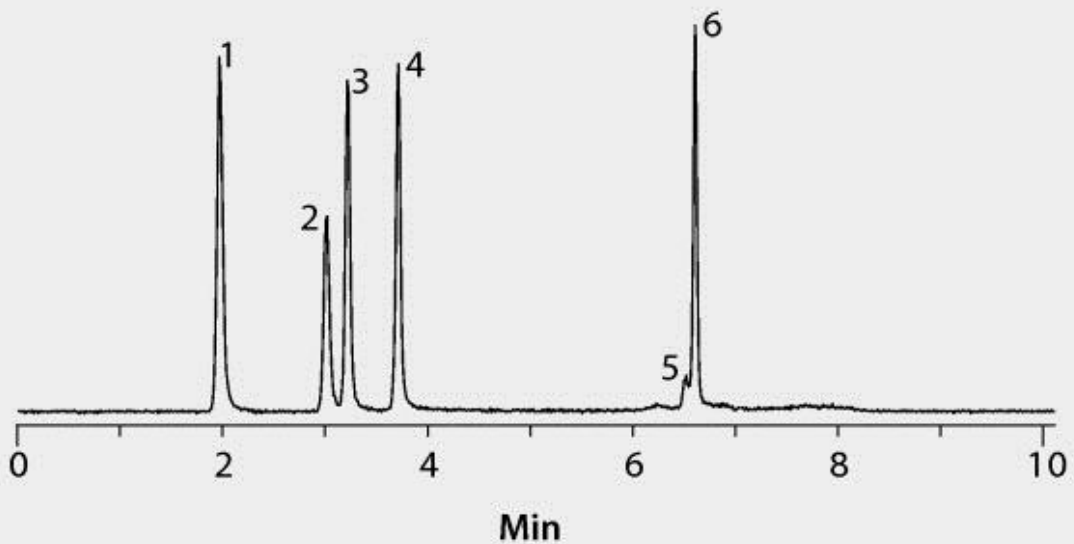

test parameters	column: Ascentis Express RP Amide, 10 cm x 2.1 mm I.D., 2.7 µm particles (53913-U)
	mobile phase: (A) 5 mM ammonium formate water; (B) 5mM ammonium formate in 95:5 (v/v) methanol:water
	gradient: held at 35% B for 0.5 min; 35 to 85% B in 5.5 min; held at 85% B for 1 min
	flow rate: 0.3 mL/min
	pressure: 320 bar
	column temp.: 35 °C
	detector: ESI(+), m/z 100-1000