


1. L-Tartaric acid
2. Lactic acid
3. Citric acid
4. Acrylic acid
5. Fumaric acid

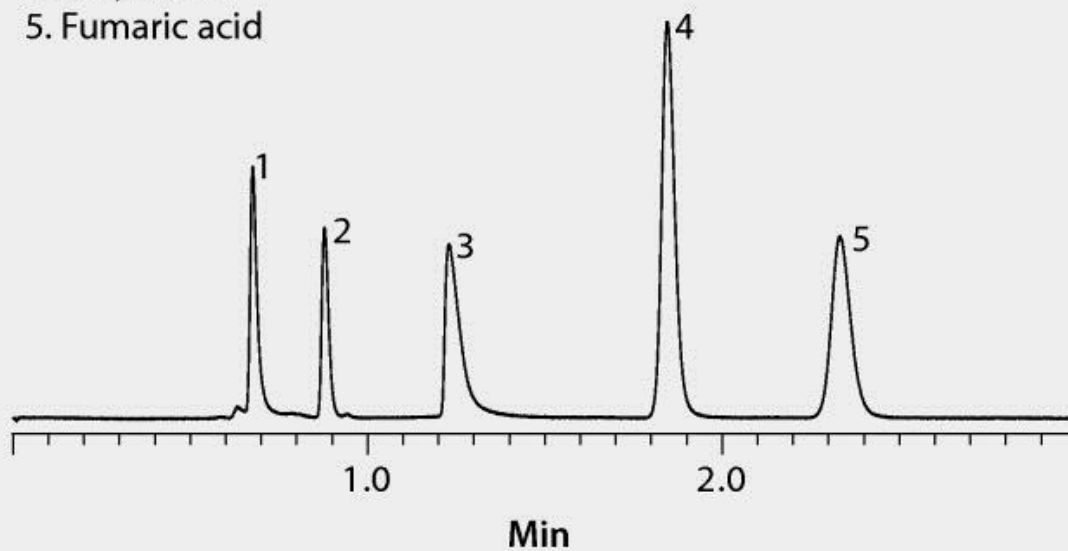

|                 |                                                                           |
|-----------------|---------------------------------------------------------------------------|
| test parameters | column: Ascentis RP-Amide, 10 cm x 2.1 mm I.D., 3 µm particles (565301-U) |
|                 | mobile phase: (A) 0.1% (v/v) TFA in (B) water                             |
|                 | flow rate: 0.3 mL/min                                                     |
|                 | column temp.: 35 °C                                                       |
|                 | detector: UV, 210 nm                                                      |
|                 | injection: 1 µL                                                           |
|                 | sample: 0.2 - 4 g/L in mobile phase                                       |