


1. PEG contaminants (m/z 520.1, 564.1, 608.2, 652.2, 696.2)
2. DAMME (m/z 604.2)

Gradient Program		
Time (min)	%A	%B
0	95	5
30	35	65

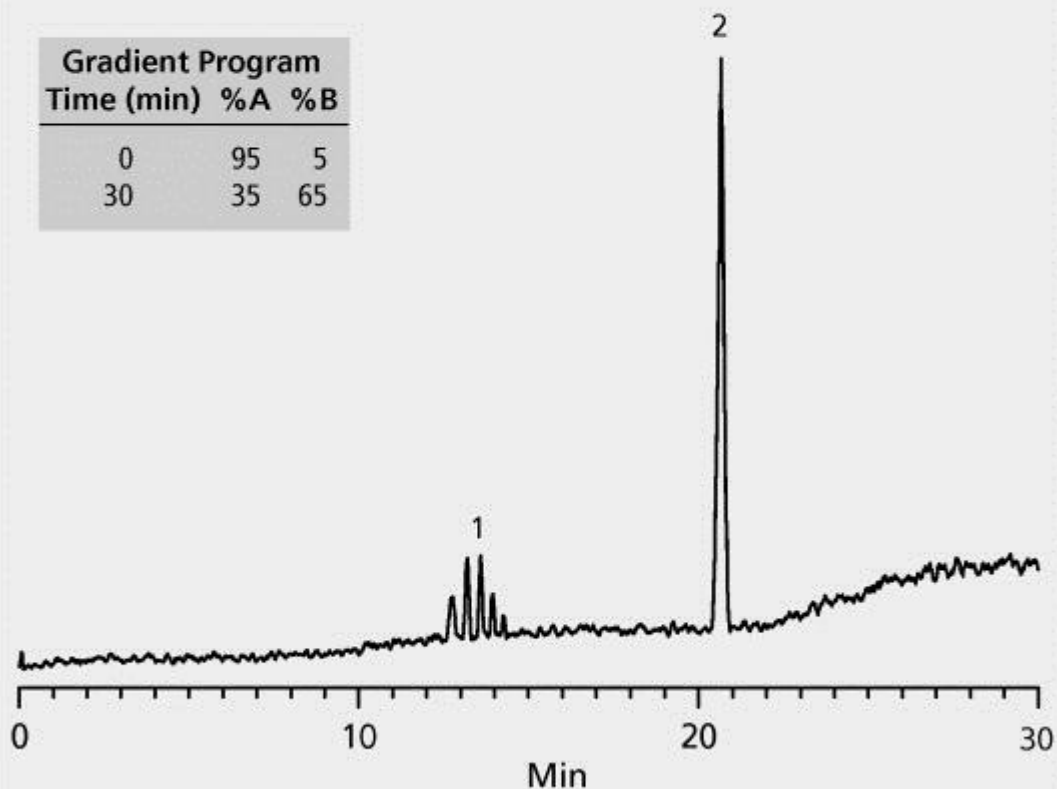

test parameters	column: Discovery HS F5, 10 cm × 2.1 mm I.D., 5 µm particles (568501-U)
	mobile phase: (A) 25 mM HCO ₂ H in water; (B) 25 mM HCO ₃ H in acetonitrile
	flow rate: 0.2 mL/min
	column temp.: 22 °C
	detector: ESI (+), full scan
	injection: 2 µL
	sample: 10 mg/L in 0.1% HCO ₂ H