




Separation of Estrogens

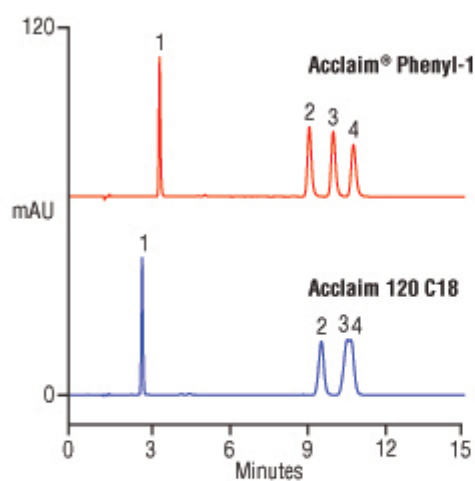


Columns: See traces
 Dimensions: 4.6 × 150 mm, 3 μm
 Mobile Phase: Acetonitrile/methanol/water
 v/v/v 15/40/45
 Temperature: 40 °C
 Flow Rate: 1.2 mL/min
 Inj. Volume: 5 μL
 Detection: UV at 220 nm
 Peaks: (100 ppm each)
 1. Estriol (E3)
 2. Estradiol (E2)
 3. Estrone (E1)
 4. Ethynylestradiol

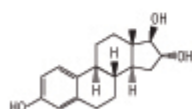


Estriol (E3)

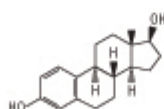

Estradiol (E2)

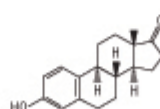

Estrone (E1)

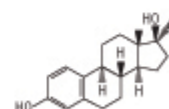

Ethynylestradiol